


Our values

INSPIRING our partner customer with quality, variety, innovation and ethics, sharing benefits and support charity .



-  +966 57 117 6000
-  info@curehealth.sa
-  Riyadh, Eastern Ring Branch Rd,
Al Ezdihar Dist 8174-3915
P.O Box 12488, KSA

Technical Specifications

Chambers	1-3
Compression mode	4
Time range	1min-24hours or continuous
Pressure range	40-150mmHg
Hold time	0-10s(Step by 1s)
Interval time	10-60s(Step by 1s)
Interface	5" touch screen
Rechargeable battery	Yes(Lithium battery :11.1V 9000mAh)
Adapter	Yes
Input voltage	AC100-240V 50-60Hz
Dimension	250(L)x184(W)x150(H)mm
Weight	2.7kg


INSPIRE COMPRESSION THERAPY DEVICE LGT-2200 DVT



LGT-2200DVT is a special Pressotherapy system that is designed to decrease the risk of DVT(Deep Venous Thrombosis) and Pulmonary Embolism(PE)before, during and after surgery

4 compression modes with reusable & disposable sleeves meet all your need

Compression Modes

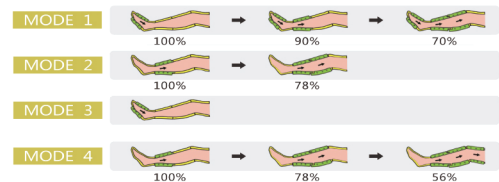

For Disposable Sleeve



*Notice: For one - leg treatment, please use single - head hose.


For Reusable Sleeve



*Notice: For two - leg treatment, please use double - head hose.

Introduction of deep vein thrombosis (DVT)

Reason

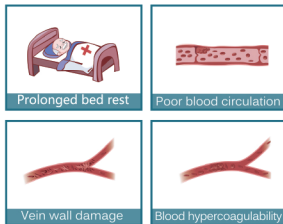

Deep vein thrombosis is one of the most common complications of bed rest, prevention is urgently needed

Features

- 4 Compression Modes
- 5" Touch Screen
- Reusable & Disposable Sleeves
- Gradient Pressure Setting
- Intermittent Pneumatic Compression
- Pressure of Every Chamber Adjustable
- Rechargeable Battery
- Powered by AC or DC



LGT-2201DVT


LGT-2201DVT is a special pressotherapy system with gradient pressure that is designed to decrease the risk of DVT (Deep Venous Thrombosis) and Pulmonary Embolism (PE) before, during and after surgery.

1 Professional Overlapping DVT Sleeve

- ⊙ Disposable Sleeve (Foot+Calf / Thigh+Calf / Foot / Calf / Thigh)
- ⊙ Reusable Sleeve



2 4 Compression Modes

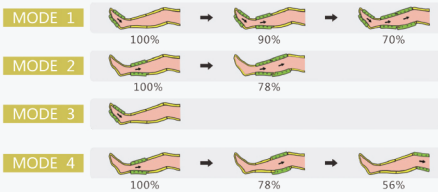


3 Gradient Pressure

Ensures the patient achieves optimal therapeutic and comfort during treatment.

4 Real-time Pressure Monitoring

The real-time pressure monitoring system can automatically fill the pressure when the pressure is lower than the setting and stop after reaching the setting.

5 Bedside / Backside Hook

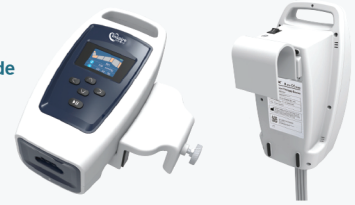

6 Flexible Treatment Time

1min-99h or continuous output

7 DC / AC Supply

8 Automatic Memory for Last Setting

9 Screen Lock Function

If you do not operate the device for about 10s, it will automatically lock to prevent misoperation.

10 Quick Connector

Securely and easily snaps into the hose



Specification

Chambers	1-3
Mode	4
Time Range	1 min - 99 hours or continuous
Pressure Range	30 - 150mmHg
Hold Time	2s (Default)
Interval Time	10s (Default)
Interface	2.2" OLED screen
Rechargeable Battery	Yes (Lithium battery: 3350mAh)
Automatic Memory	Yes
Dimension	270(L) x 148(W) x 145(H)mm
Weight	2.2kg

⊕ Feature

- ⊗ Made of premium fabric with more lightweight, breathable and comfortable
- ⊗ Compatible with both Inspire™ LGT-2200DVT and LGT-2201DVT
- ⊗ With a wide range of types for different application
- ⊗ One-Three (1-3) chambered sleeve with intermittent and sequential (gradient) inflation

DVT Sleeve Size Details

Order No.	Sleeve	Size	Dimension	Proximal Circumference
220092M	Thigh Length Sleeve	M	695 (L) x 557(W) mm	<69.5cm(28")
220092L		L	914 (L) x 649(W) mm	<91.4cm(36")
220093M	Calf and Foot Sleeve	M	587.6(L)x670.4(W) mm	<58.7cm(23")
220093L		L	under pending	under pending
220096	Thigh Sleeve	S	750 (L) x 250(W) mm	65-75cm (26-30")
220096M		M	840 (L) x 275(W) mm	74-84cm (30-34")
220096L		L	960 (L) x 312(W) mm	86-96cm (35-39")
220055	Calf Sleeve	S	620 (L) x 366(W) mm	50-62cm (20-24")
220055 M		M	717.1 (L) x 378.5(W) mm	60-71cm (24-28")
220055 L		L	812.6 (L) x 416.7(W) mm	70-81cm (28-32")
220095	Foot Sleeve	S	443.7 (L) x 152(W) mm	/
220095M		M	556 (L) x 167.4(W) mm	/
220095L		L	693 (L) x 184.3(W) mm	/



Thigh Length Sleeve



Calf and Foot Sleeve



Calf Sleeve



Foot Sleeve

Specifications

- ⌚ Channel: 4
- ⌚ Chambers: 4/8
- ⌚ Compression modes: 8
- ⌚ Time range: 1-99 min
- ⌚ Pressure range: 37-187 mmHg
- ⌚ Interface: 5.7" Touch screen
- ⌚ Weight: 10 kg
- ⌚ Dimension: 366mm(W) x 393mm(L) x 194mm(H)



Features

- ⌚ User-friendly interface, just a few steps to finish your setting
- ⌚ Gentle inflation and deflation
- ⌚ Durable & easy-use sleeves with overlapping chambers
- ⌚ Pressure unit switchable between mmHg and kPa
- ⌚ 8 compression modes, 4 output channels
- ⌚ Real-time pressure, and inflation part monitoring
- ⌚ 2 pcs of 8-chamber cuffs can be connected with one device at the same time

Application

- ⌚ Lymph drainage
- ⌚ Blood circulation
- ⌚ Detoxification
- ⌚ Cellulite treatment
- ⌚ Body shaping in legs, abdomen, buttocks and arms
- ⌚ Strengthening of loose skin
- ⌚ Relaxation of stress
- ⌚ Relieve fatigue

Mechanism & Functions

The multi-lymphatic drainage system LGT-2200WM is used to drain the dissolved or melted fat cells generated after various obesity equipment by improving lymph glands or blood vessels that are constricted or damaged due to cellulite or excessive fat cells on the lymphatic drainage step, one of the general obesity lipolysis steps.



■ Deep Vein Thrombosis (DVT)